


New Product Update

MARCH/APRIL 2022

INNOMED

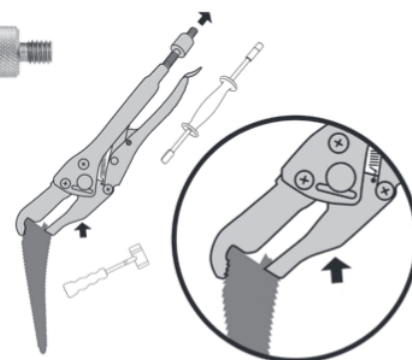

PRODUCT NO'S:

3976-00	[Broach Extraction OrthoVise Set with Small Slaphammer]
Set Includes / Available Individually:	
3976-01	[Broach Extraction OrthoVise Only] Overall Length: 9" (22,9 cm)
3955	[Slap Hammer for Small OrthoVise] Overall Length: 8.75" (22,2 cm)
3985-03	[Threaded Adapting Screw-Small]

Broach Extraction OrthoVise™

Designed by Joel Matta, MD

Designed for hip broach extraction when the broach post is broken or there is a failure of the broach handle

Hand/Finger Positioner

Designed by Emad Aboujaoude, MS, MPAS, PA-C

Designed to help provide surgical positioning during fluoroscopy and fixation by isolating the operative digit while retracting the unaffected digits

Radiolucent positioner can be steam or gas sterilized.

PRODUCT NO:

1134
Overall Length: 5.75" (14,6 cm)
Handle Width: 4.25" (10,8 cm)
Blade Width: 1.6" (4 cm)



Uses include but not limited to:

- ▶ Intramedullary Metacarpal Screw
- ▶ Phalanges CRPP
- ▶ Digit Amputation
- ▶ Digit Mass Excision
- ▶ Finger Joint Fusion



Prototype Shown

Modified Angled Hohmann Retractor with Long Handle and Short Tip

Designed by R. Michael Meneghini, MD


Longer handle to help provide safe patella retraction with excellent ergonomics, and useful in other orthopedic procedures



PRODUCT NO:

7119
Overall Length: 11.5" (29,2 cm)
Handle Length: 9.25" (23,5 cm)
Blade Depth: 5" (12,7 cm)
Blade Width: .75" (19 mm)
Prong Depth: 1" (2,54 cm)

Profile View

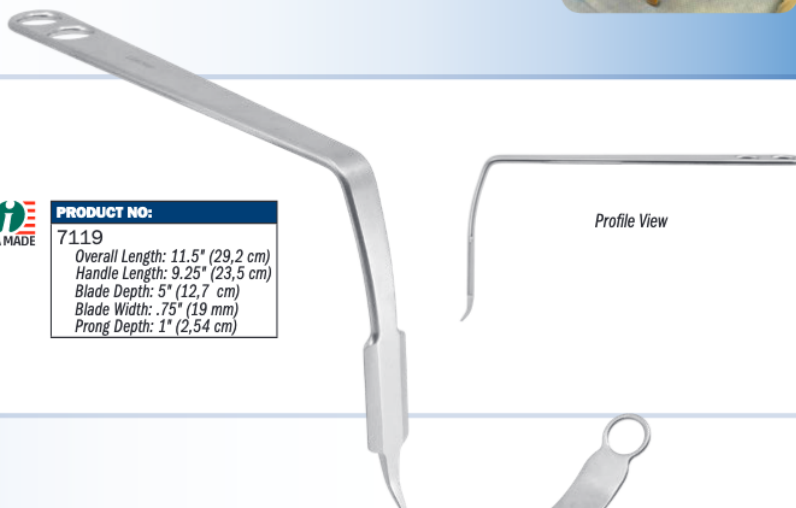

Modified Wide Hohmann Retractor with Taylor Tip

Designed by Jeffrey P. Beckenbaugh, DO

Anterior and posterior acetabular retractors for all approaches, including the direct anterior approach, featuring a hammer platform for insertion with a mallet

Used as a calcar and posterior femoral retractor for the posterior approach, and an anterior femoral elevator for the direct anterior approach.

PRODUCT NO'S:

3012-00	[Set - Left & Right]
Also available individually:	
3012-L	[Left] Overall Length: 14.75" (37,5 cm) Blade Depth: 5" (12,7 cm) Blade Width: 45 mm
3012-R	[Right] Overall Length: 14.75" (37,5 cm) Blade Depth: 5" (12,7 cm) Blade Width: 45 mm



McMaster Abductor Retractor

Designed by William D. McMaster, MD

Designed to help with proximal femur exposure helping to protect the abductors – gluteus medius and minimus – during posterior approach THA

The ergonomic design allows any application where soft tissue retraction is needed for exposure similar to a doyen or hibbs-type retractor.

PRODUCT NO:

6385


Overall Length: 8.75" (22.2 cm)

Handle Length: 4.625" (11.7 cm)

Blade Depth: 4.125" (10.8 cm)

Blade Width: 1.5" (3.8 cm)

Blade Curve Diameter: 1.75" (4.4 cm)

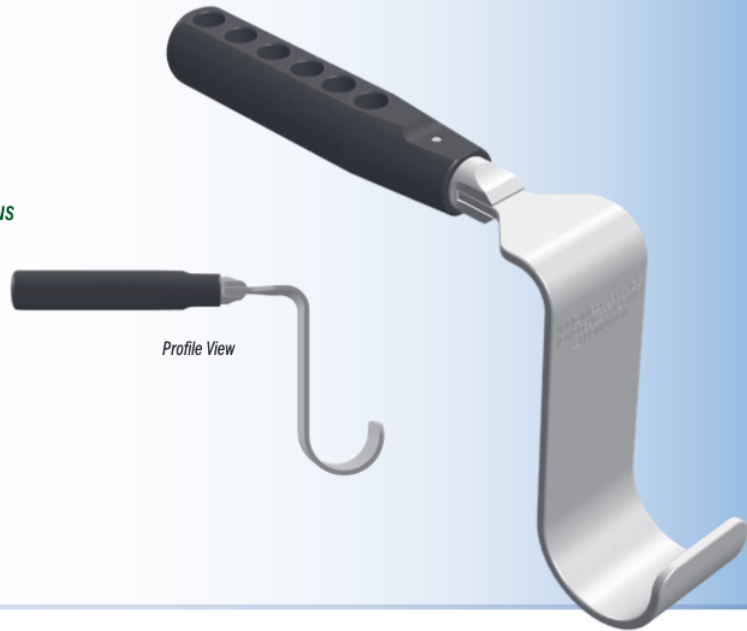

Profile View



Retractor Clip for Smoke Evacuation Tube

Designed by James Saucedo, MD

Repositionable stainless steel fastener designed to clip onto a retractor to help control the location of a smoke evacuation tube

Allows for use on a 1/8" thick material with allowance for a "spring" fit.

PRODUCT NO:

5466



Screw Extractor with Speed Lock

Designed by Khaled Sarraf, MD & Konstantinos Doudoulakis, MD

Universal extractor designed to accommodate a large range of screws and screw heads

Can also be used to help with removal of other devices that may require a twisting universal locking gripper.

PRODUCT NO:

2021

Overall Length: 9.25" (23.5 cm)

Jaw Width: 11 mm

Jaw Length: 5 cm



INNOMED

© 2022 Innomed, Inc.

Innomed, Inc
103 Estus Drive
Savannah, GA 31404
USA

Tel 912.236.0000
Fax 912.236.7766
www.innomed.net
info@innomed.net

Innomed-Europe LLC

Alte Steinhäuserstrasse 19
CH-6330 Cham, Switzerland
Tel 0041 (0) 41 740 67 74
Fax 0041 (0) 41 740 67 71

www.innomed-europe.com
info@innomed-europe.com
contact.france@innomed-europe.com

Innomed-Europe GmbH

c/o Emons Logistik GmbH
In Rammelswiesen 9 D-78056
Villingen-Schwenningen
Deutschland

Tel 0049 (0) 7720 46110 60
Fax 0049 (0) 7720 46110 61
www.innomed-europe.com
info@innomed-europe.com